Showlon Portable Irrigator

User Manual

Please read the manual carefully before using


About Oral Irrigator

Oral Irrigator is an auxiliary tool for cleaning oral cavity. It can clean the hidden positions such as the side gap of teeth and the gap between gum and tooth by pulse water washing the oral cavity, at the same time, it can massage and stimulate the gum, so as to improve the health of teeth. S3 adopts portable storage design, which is light and compact. It is easy to carry at home, office and travel. It has three modes: soft, pulse and strong. It can meet different personal needs and adapt easily. S3 can be used as a companion to take care of your teeth.

Notes

- Please choose the right mode to avoid discomfort of teeth and gums
- Please fill purified water as much as possible. Don't fill iodine, alcohol, salt water, medicine or other essential oil insoluble in water into the water tank.
- Don't fill the tank with water above 40℃
- Don't use the product to clean nose, ear or other vulnerable and sensitive body parts than oral cavity, since water pressure may cause injuries.
- For patients with periodontal disease or who have had oral surgery in the last two months, please consult your doctor before using.
- Try to avoid sharing this product with others
- Please keep the charging port dry
- When storing, please be careful to clamp hands at the joint between the water tank and the shell

Product Layout

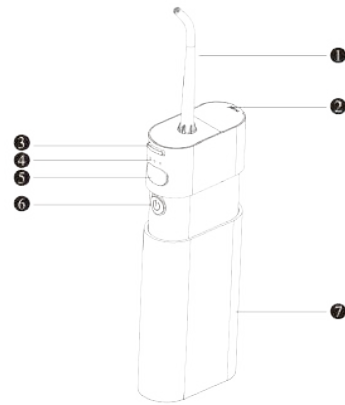- In The Box:
S3*1
Nozzle*1
User Manual*1
- Nozzle
 - Water Inlet
 - Pull out and Insert Button of Nozzle
 - Indicator Light
 - Charging Port
 - On/Off Button
 - Water Tank

How to use

1. Install the nozzle
2. Open the water tank and the water inlet to fill water
3. Put the nozzle into mouth, press the on/off button to use
4. Long press the on/off button to switch the modes
5. Pull down the nozzle, open the water inlet cover and store the nozzle
6. Pull out the tank and clean it

Modes Description

Modes	Mode Switch	Indicator Light	Target User	Water Pressure
Soft Mode	Tap the on/off button	Blue	The water pressure is gentler and suitable for sensitive teeth and first-time users	87-100PSI
Pulse Mode	Long press the on/off button 0.5s	Blue	Variable-frequency pulse water pressure, strong and massage gums	70-85PSI
strong Mode	Long press the on/off button 0.5s again	Blue	Strong water pressure, high efficient cleaning, meet the daily cleaning for most people	≥110PSI
Indicator Light	1. When power is less than 10%, it turns blue 2. When charging, it turns flashing blue 3. When full charged, it turns blue, and turn off after 2 minutes			

Tips:

Gums may bleed or itch after first time use, it is recommended to use pulse mode or soft mode, and then use strong mode after you get used to it.

Product Specification

Product Name	Mini Portable Oral Irrigator	Model	S3
Nozzle Diameter	6.6mm	Pulse Frequency	1600 times/minute
Water Pressure	70-110PSI	Water Tank Capacity	110ml
Battery Capacity	1100mAh	Noise	75dB
Input	5V=1A	Rated Power	6W
Charging Port	Type-C	Charging Time	About 4h
Waterproof	IPX7	Modes	3
Product Size	130*63*30mm	Net Weight	175g

Basic Maintenance

1. Nozzle
 - Press the pull/insert button of the nozzle and remove it, clean it with water and then wipe it with a rag
 - Dentists suggest replacing nozzles every six months
 - If the top of nozzle is deformed, please replace it immediately to avoid product breakdown
2. Body
 - Remove the body dirt with a cloth dipped in detergent, flush it with clean water and wipe it clean with a dry cloth
 - Don't use diluent gasoline alcohol and other solvents to wipe the product
 - Don't clean the product with water above 50℃
3. Tank
 - Pull out the tank regularly, remove the inside dirt with a brush and wipe it clean with a dry cloth
 - Dry out the tank if not using for a long time

Troubleshooting

Trouble	Probable Reasons	Solutions
Not working	Battery low	Charge the product
Indicator light abnormal/ not on	Bad connection to the charging cable	Unplug the cable and plug again
	Faulty charging circuit	Contact the seller
Short use time after charging	Insufficient charging time	Charge it more than 4hours
	Faulty circuit	Contact the seller
Low water pressure	Nozzle block or damage	Replace the nozzle

Warranty and Service

In order to protect your rights and interests, please read the warranty regulations carefully, you will enjoy the after-sales service what we provide. We promise:

- Within 12 months from the date of purchase, if the product has quality problems and meets the scope of free maintenance, please keep the product and accessories complete. We will provide free maintenance
- Please refer to the purchase platform for the specific time limit and process of return and exchange. You can also contact the seller directly for after-sales service

During the warranty period, we will not provide after-sales service under any of the following circumstances

- Faults caused by human, including those caused by improper use, maintenance, storage or operation not in accordance with the instructions, as well as those caused by self disassembly and assembly of products
- Failure caused by force majeure (such as fire, earthquake, flood, etc.)
- Unable to provide purchase certificate or modify the purchase certificate without permission

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